

From Eye to Insight

Leica
MICROSYSTEMS

LEICA M620 F20

Easy-to-use, high quality surgical microscope
for Ophthalmology



LEICA M620 F20



The Leica M620 surgical microscope is a perfect balance of form and function. Precision optics combined with the finest mechanics make the Leica M620 an outstanding choice. Swiss design, quality and precision – the Leica M620 meets the needs of modern day ophthalmic surgery.

High-quality and reliable performance

The Leica M620 features OptiChrome optics from Leica Microsystems and a direct halogen illumination system. These features deliver the quality and reliability of a premium ophthalmic microscope at an affordable price.

Compact, light and easy to move

The Leica F20 floorstand is compact and moves easily. With the high-quality precision bearings and a long swingarm, the stand is highly maneuverable and easy to set up for surgery.

Leica Design by Christophe Apothéoz

BRILLIANCE IN OPTICS AND LIGHT



Leica OptiChrome for higher contrast and crisper, sharper images in natural color.

Leica OptiChrome Optics

Leica Microsystems' crisp, sharp OptiChrome optics give the surgeon natural color, outstanding depth of focus, and high contrast for enhanced detail recognition. The high degree of light transmission with OptiChrome technology allows the surgeon to use lower degrees of illumination and thus supports patient safety.

Direct Halogen Illumination

Leica Microsystems' exclusive Direct Halogen Illumination system ensures that the surgeon will have homogeneous illumination to see crisp images even at low light levels. Direct illumination is a feature that consistently places Leica ophthalmic microscopes above the systems with fiber optics illumination.



Under sterile conditions, the jalousie provides outstanding Red Reflex.

Outstanding Red Reflex

Based on more than 25 years of experience with direct halogen illumination, Leica Microsystems developed an illumination concept that delivers bright, stable Red Reflex. Even with difficult starting positions or heavy pigmentation of the patient's eye, the Red Reflex has outstanding intensity.

Individually optimized Red Reflex

The jalousie, which can be operated in sterile conditions, allows perfect illumination for any surgical situation with optimum contrast and depth of field. It helps to further enhance and optimize the Red Reflex.

Proven, Direct Halogen Illumination
for crystal-clear images



INTUITIVE CONTROL AND INCREASED SAFETY



Intuitive controls for light intensity, focus, zoom, and XY speed.

Straightforward operation

Leica Microsystems' touchscreen, built into the floorstand for convenient, easy access, offers intuitive control of all Leica M620 functions. The surgeon can see the microscope's functions at a glance on the screen. Individual user settings can be programmed for each surgeon to increase efficiency.

Auto Reset and the microscope is ready-to-use

When the surgeon pushes the swingarm up after each procedure, all microscope functions automatically reset and the microscope is ready for the next case.

Antimicrobial coating for added safety

Leica surgical microscopes are designed to support your clinical needs including your hygiene protocols. Our microscopes are therefore coated with a paint which is designed to provide an antimicrobial effect on surfaces.



FLEXIBILITY FOR INDIVIDUAL NEEDS



Leica Rotatable Beamsplitter

Leica M620 with Video Adapter at the rear and rotatable assistant microscope tube.

Leica Microsystems accessories, including video cameras, adapt seamlessly with the Leica M620 series microscopes.

Going temporal? Swing around...

For those surgical procedures where an assistant is needed, especially with the growing number of temporal approach cataract surgeries, Leica Microsystems proudly presents the world's first two-beam path solution – the Leica Rotatable Beamsplitter. The easy side-to-side quick change of the assistant observer optics enables time-saving between cases, increases efficiency in the operating room and minimizes potential damages to the binoculars. Also the interface for the Leica Video Adapter is conveniently placed at the rear for maximum flexibility.



Standard Beamsplitter

Leica M620 with Video Adapter and assistant microscope tube. The stereo assistant tube can be attached to either side of the beamsplitter.

Wide variety of binoculars provide comfort

Long working hours under the wrong ergonomic conditions can cause fatigue and discomfort. For this reason Leica Microsystems has a wide selection of binoculars and objective lenses to choose from to provide comfort for the surgeon.



Leica M620 with new 5°–25° binocular tube and compact handles.



Leica Microsystems' range of binocular ergonomic solutions to meet individual comfort needs.

FLEXIBILITY FOR INDIVIDUAL NEEDS



Unique Leica Zoom Video Adapter

The Leica Zoom Video Adapter with its 3× optical zoom and fine focus allow individual magnification adjustment of the video image separate and independent of the microscope zoom. A sterile fine focus adjustment on the video adapter fine tunes images to perfection.

Integrated Leica 2D Video System

The Leica 2D video module can be easily integrated with the beam path of the Leica M620 F20 and offers good ergonomic conditions for the surgeon due to its compact design.

Retinal accessories

For application in posterior segment surgery the Leica M620 can carry the most common accessories like inverters, Oculus BIOM system or laser-filters.

The wide selection of foot switches fulfills any ergonomic need. All functions can be individually programmed according to the preferences of each surgeon.

Leica Wireless Footswitch

Using ISM bandwidth technology, it offers maximum mobility for a fast and easy switch between left and right eye procedures.



TECHNICAL SPECIFICATIONS

OPTICS LEICA M620

| | |
|-----------------------|---|
| Magnification changer | 6:1 zoom, motorized, adjustable speed |
| Focusing | Motorized, 50 mm, adjustable speed, with automatic reset |
| Objective | APO Leica WD = 175 mm, 200 mm Achromatic Leica f = 175 mm, 200 mm and 225 mm |
| Eyepieces | Widefield eyepieces for those wearing glasses 8.33×, 10× and 12.5×; Diopter adjustment +/- 5 with adjustable eyecup |
| Illumination | Homogenous, coaxial illumination; Integrated Red Reflex illumination, light intensity adjustable via control unit and/or footswitch jalousie for continuous, situation-dependent reduction of ambient illumination via sterile controls; Quick-change lamp mount with 2 halogen lamps, 12V/50W filter changer with UV projective filter; 2 free openings for selectable filters diameter 32 mm |
| Inclining mechanism | +/- 5° manual fine adjustment, sterile operation |
| XY-unit | Positioning range 50 mm in X and Y direction; Adjustable speed, reset button, automatic reset at power on |
| Control unit | Multifunctional display for brightness, magnification, reset and user selection, adjusting the speed of drives, individual user programming, Auto Reset |
| Hand/Footswitch | 12-function footswitch with controls arranged longitudinally or transversally, wired or wireless Optional: 6- and 16-function wired footswitch, 12- and 16-function wireless footswitch, 12-function handswitch |

FLOOR STAND LEICA F20

| | |
|--------------------------|--|
| Type | F20 Compact floorstand |
| Base | Compact base with 4 movable rollers, 4 foot brakes |
| Base dimension | 608 × 608 mm |
| Maximum load | 6.5 kg accessories on the microscope |
| Maximum reach | 1444 mm |
| Minimum transport height | 1902 mm |
| Weight approx. | 210 kg (minimum equipment) |


ACCESSORIES

| | |
|------------------------|--|
| Assistant's attachment | Stereo co-viewer, Monocular co-viewer |
| Beam splitter | 50% / 50%, 70% / 30% Leica Rotatable Beamsplitter |
| TV/photo | Leica 2D video system Leica Zoom Video Adapter (ZVA), 3:1 zoom, 35 mm to 100 mm focal length, c-mount, with fine focus Leica Manual Video Adapter (MVA), 55 mm, 70 mm focal length, c-mount, with fine focus Photo/TV dual attachment f = 60/85/107 mm for TV, f = 250/350 mm for 35 mm camera Photo attachment f = 250/350 mm |
| Wide-angle observation | BIOM, EIBOS |
| Inverters | AVI, SDI, OIVSL, ROLS |
| Laser | Can be fitted, adapters available from laser suppliers |
| Asepsis | Sterilizable components for all controls |

ELECTRICAL DATA

| | |
|------------------|--|
| Supply | 100–240 VAC (+/- 10%), 50/60 Hz / 250 VA |
| Protection class | Class 1 |
| Safety type | Type B |

-

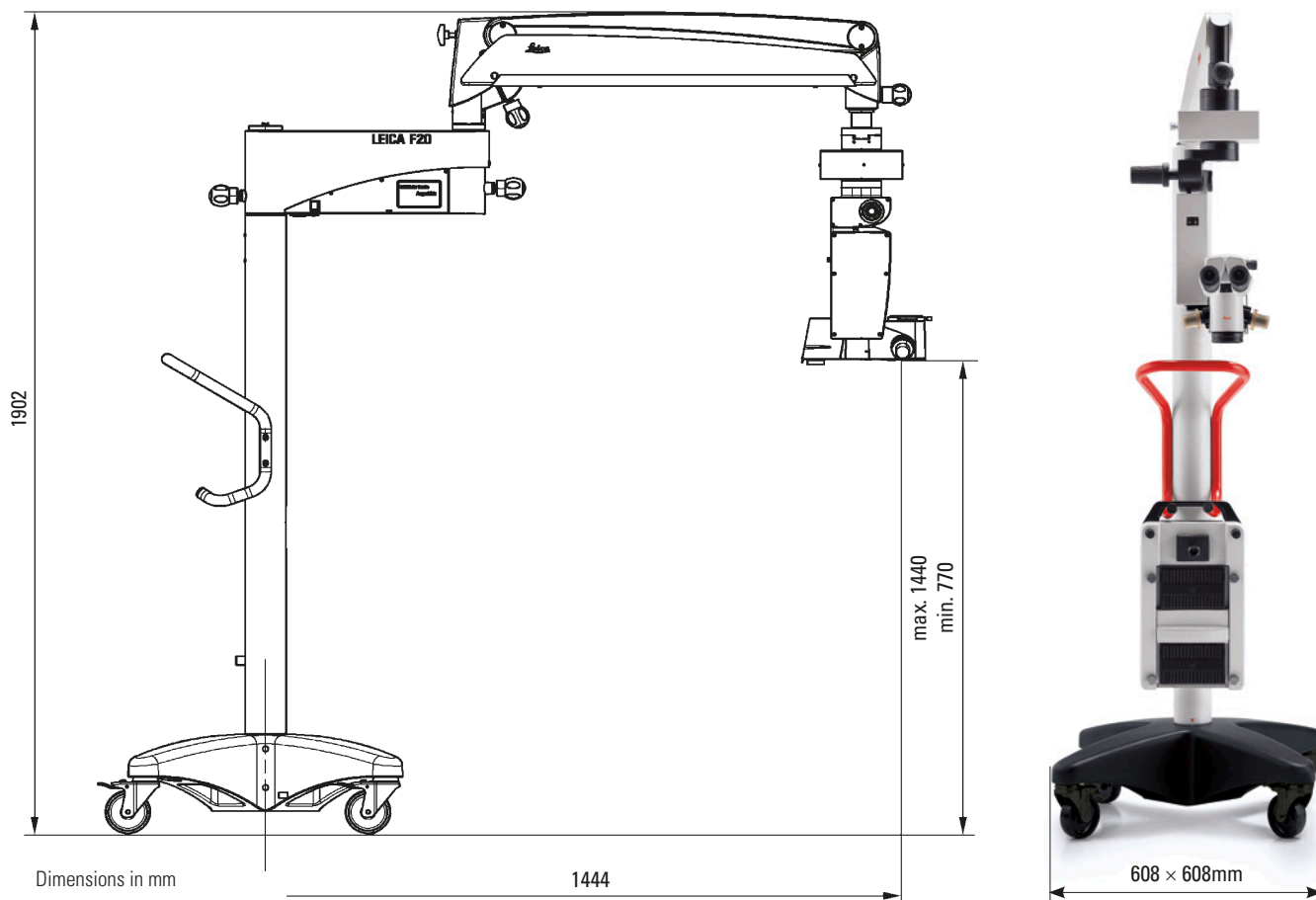

Leica Microsystems (Schweiz) AG
Max Schmidheiny-Strasse 201
9435 Heerbrugg, Switzerland



Class I surgical microscope M620 F20 incl. accessories.

Not all products or services are approved or offered in every market and approved labeling and instructions may vary between countries. Please contact your local Leica representative for details.

Leica Microsystems (Schweiz) AG · Max Schmidheiny Strasse 201 · CH-9435 Heerbrugg

T +41 71 726 3333 · F +41 71 726 3399

www.leica-microsystems.com

CONNECT
WITH US!

